


# EtchMaster® Blows Away the Competition

## Ultimate Powder Blasting Tips for Intra-Oral Procedures Prophylaxis & Abrasion

### Perfect for:

- surface preparation prior to sealants
- composite & porcelain repairs
- crowns, veneers, brackets, bonding
- **PLUS** localized stain removal & surface polishing.

### No Mess – Costs Less

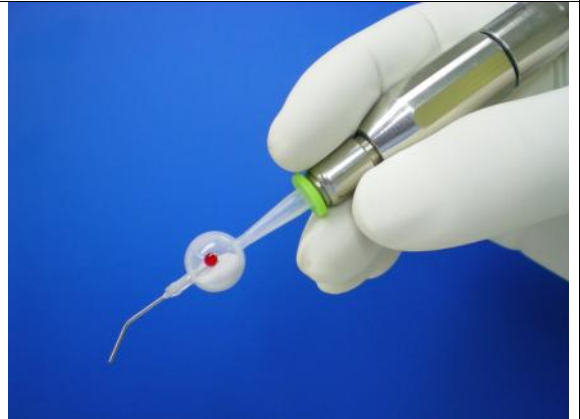


- Unidose • Pre-Filled • No Clean-up • No Equipment • No Overspray • Never Clogs
- Simple Adaptors to Chairside Air • 10 Free Tips with every Adaptor

### Comparison to other systems

Clinically effective and mess-free using the same powders, air pressures, and nozzle sizes of traditional non-disposable air abrasive equipment.

### Advantages over other systems

- No Overspray • No Clogging • No Equipment/Maintenance • No Refilling • No Nozzle Replacement • No Cleanup

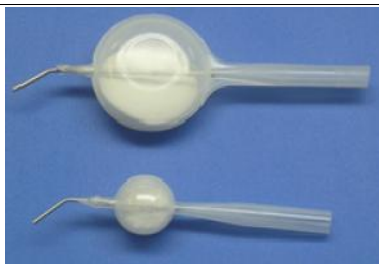

### Award winning NO MESS™ pre-filled disposable abrasive tips

improve cementation bond strength to tooth surfaces, prosthetics, and sealants for intra-oral General, Cosmetic, Paediatric, and Orthodontic treatments; Mild and soft abrasive EtchMaster® tips provide localized stain removal and prophylaxis capabilities. Simple and inexpensive EtchMaster® adaptors are a cost effective replacement to traditional messy

micro air abrasive etcher systems.

The EtchMaster® tips offer you a Carefree™ and No Mess™ Air Abrasion experience by solving the deficiencies of traditional air abrasive equipment, saving you Time & Money:

### How do I connect?

Adaptor screws directly to the standard chair-side handpiece connector – universal compatibility (3 – 6 hole)

Air flow is controlled via the standard chair-side foot pedal

Air flow is regulated to optimum range of 35 – 45 psi.



# EtchMaster® Adaptor Selection



## Pressure Regulated 4-Hole Adapter £490.00

PR Adaptor – self regulates to optimal air pressure regardless of chair airline setting.

(Note: water delivery port is sealed off)

Adaptor screws directly to the standard chair-side handpiece connector – universal compatibility (3-6 hole). Air flow is controlled via the standard chair-side foot pedal. Air flow is regulated to optimum range of 35-45 psi.




## EtchMaster® Kavo Style Adapter Starter System £399.99

Adapter mounts directly to Kavo handpiece connector

Air flow is controlled via the standard chairside foot pedal

A Pencil Style handpiece provides easy reach

(Note: water delivery port is sealed off)

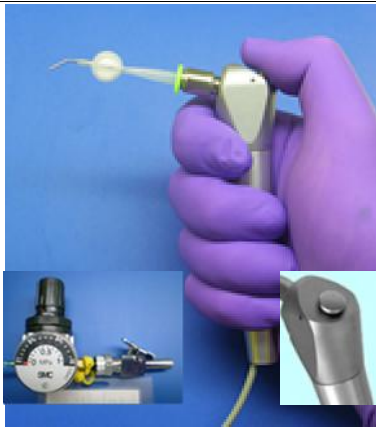

## Quick Disconnect Adaptor £399.99

This is a direct replacement for MicroEtcher type systems.

Adapter plugs into a Quick Disconnect (QD) type connection. Integral to the Quick Disconnect adapter is a miniature pressure regulator. The pressure regulator is pre-set to 40psi eliminating any patient sensitivity.

Air flow is controlled via a syringe type finger actuation

(Note: adapter is pressure limited so it may be connected to a high pressure air source without over-etching)



## Pre-filled colour-coded disposable abrasive tips:


### Two Tip Sizes

- Small size tips cover 4 tooth surfaces £2.19 inc. VAT
- Large size tips cover 16 tooth surfaces £3.35 inc. VAT

### Small Tips Usage:

Sealants for 4 molars, all cosmetic bonding, bracket repair, composite & porcelain repairs, and localized polishing

**Large Tips Usage:** Full mouth sealants, crown prep, orthodontic brackets, stain removal and polishing.

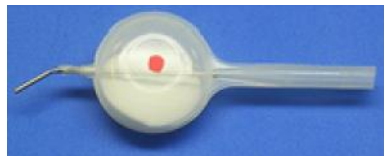

**Available Powders:** 27 micron & 50 micron Aluminium Oxide, Glass Beads, Aluminium Oxide & Sodium Bi-Carbonate, Calcium Carbonate, Sodium Bi-Carbonate.

Colour-coded tips

To order call DPS on: 01438 820550

Email: [info@d-p-s.uk.com](mailto:info@d-p-s.uk.com) [www.d-p-s.uk.com](http://www.d-p-s.uk.com)

